

Urgent Field Safety Notice
Deep Brain Stimulation Therapy
Reported Events Related to DBS Tunneling Procedure
Labeling Update

May 2016

Medtronic reference: FA712

Dear Healthcare Professional,

The purpose of this letter is to notify you of information that will be added in the future to Medtronic's Deep Brain Stimulation (DBS) labeling, in connection with the tunneling procedure related to DBS Extensions Models 7483, 37086 and 37087. These labeling updates result from Medtronic's ongoing monitoring of reported events and published literature. The labeling updates further clarify potential risks which have been reported with the DBS implant procedure. Medtronic is sharing this information to help you with the DBS implant procedure, patient selection, informed consent, and post-implant follow up for patients treated with DBS Therapy.

Background:

Medtronic's post marketing surveillance for DBS therapy has identified the following serious events or injuries associated with tunneling that are not presented in the current DBS labeling. The following events are associated with the DBS implant and tunneling the extension from the lead to the INS:

- Cases of spinal accessory nerve injury
- Severe bleeding or vascular injury
- An extension inadvertently implanted through the ribs and later found to be close to the heart requiring revision
- Tunneling too superficially resulting in an extension that exited and then re-entered the neck which wasn't discovered until the surgical drape was removed

These events occurred in an estimated 0.008% of tunneling procedures. Note: The specific model of tunneling tool associated with the patient injuries was not reported in many of these events.


Recommendations:

During the DBS implant, use caution while tunneling the extension from the lead to the INS to avoid tunneling too deeply or superficially because serious injury may occur. Complications or effects related to the tunneling procedure may include injury to nerve tissue (such as the spinal accessory nerve), vascular injury that may result in prolonged hospitalization, and tunneling through unintended anatomy (such as in between the ribs and entering the thoracic cavity).

The Competent Authority of your country has been notified of this action.

Thank you for your continued support. We apologise for any inconvenience this may cause. We remain committed to continuing to improve our product performance and services to enable you to manage your patients in a safe and effective manner. In case you have any questions regarding this notification please contact your Medtronic representative Nynne Stewart at 2672 6815 or at nynne.stewart@medtronic.com.

Sincerely,

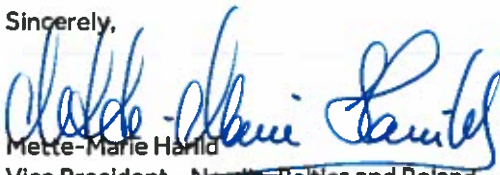

Mette-Marie Harild
Vice President – Nordic, Baltics and Poland