

Instruction clean probe assembly cobas t 711

Difference between old and new procedure

cobas t 711		
Step	Old procedure	New procedure
1	Open Main Cover	Open Main Cover
2	Remove drip trays	Remove drip trays
3	Clean probes	Clean probes
4	Clean drip trays	Clean drip trays
5	Reattach drip trays	Reattach drip trays
6	Move reagent transfer arm	Move reagent transfer arm
7	Remove tube holders	Moisten the cleaning swab with disinfectant, and clean the inside of the tube holder
8	Clean tube holder	
9	Reattach tube holder	

⚠ WARNING

Contact with sample probes


Contact with sample probes may result in injury and infection.

- ▶ Avoid touching the end of the sample probes
- ▶ Avoid contacting the sample probes during maintenance
- ▶ Wear appropriate personal protective equipment.
Take extra care when working with lab gloves. They can easily be pierced or cut, leading to infection.

Parts Required

SWABS LINT-FREE Mat#08288267001 (Set of 50), Chemtronics® Foamtips #120 Swabs

Note: This material is also part of the startup kit of the **cobas t 711**

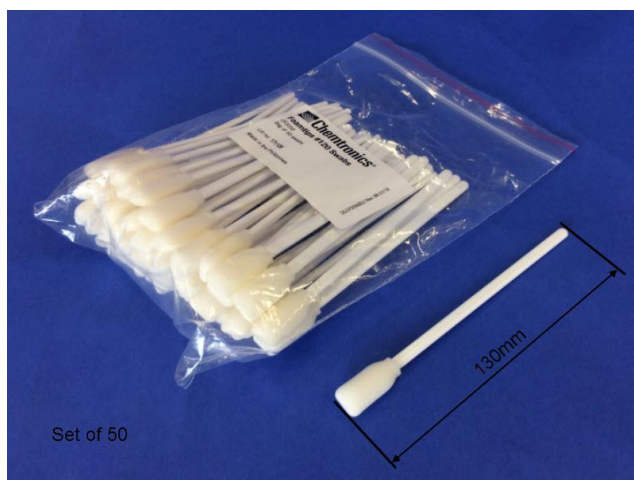

Instruction clean probe assembly cobas t 711

Instructions Maintenance Action 'Clean probe assembly'

As usual start 'Clean probe assembly' maintenance action in the software, then follow steps as described below also in the software wizard.


1. Open main cover and click *Next Step*



Open the main cover.

WARNING: Risk of infection due to contact with contaminated surfaces. Wear lab gloves.

2. Remove drip tray and click *Next Step*

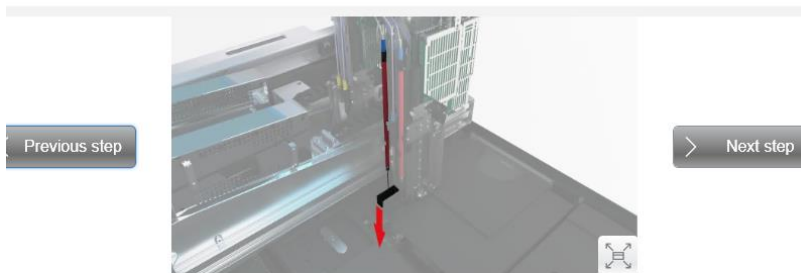


While ensuring that the probes are not bent, carefully remove the drip trays in a downward motion beginning with R1.

WARNING: Risk of personal injury and infection. Use caution around the reagent probes. Wear lab gloves.


Instruction clean probe assembly cobas t 711

3. Clean probes and click *Next Step*

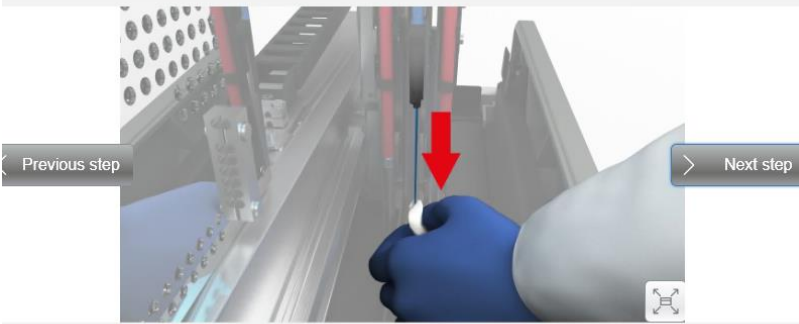

While ensuring that the probes are not bent, carefully clean the probes from top to bottom with a lint-free cloth moistened with disinfectant.


4. Clean drip trays and click *Next Step*

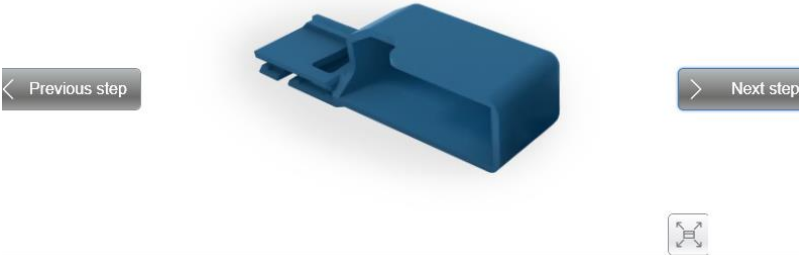


Clean the inside and outside of the drip trays.

Instruction clean probe assembly cobas t 711

5. Reattach drip trays and click *Next Step*

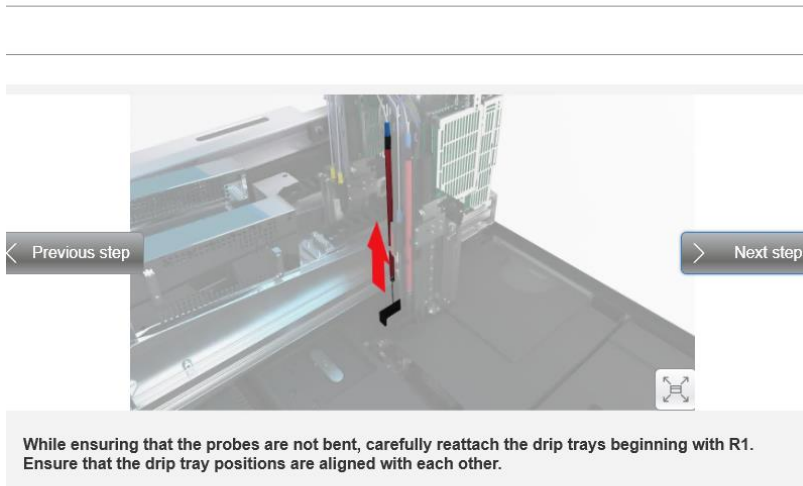


6. Move reagent transfer arm and click *Next Step*

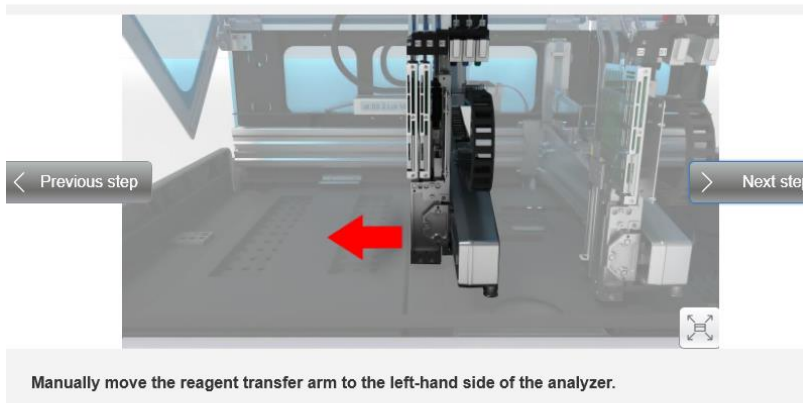


Important: Do not continue to remove tube holder, continue now with step 7 below.

Instruction clean probe assembly cobas t 711

7. Moisten the cleaning swab with disinfectant, and clean the inside of the tube holder.

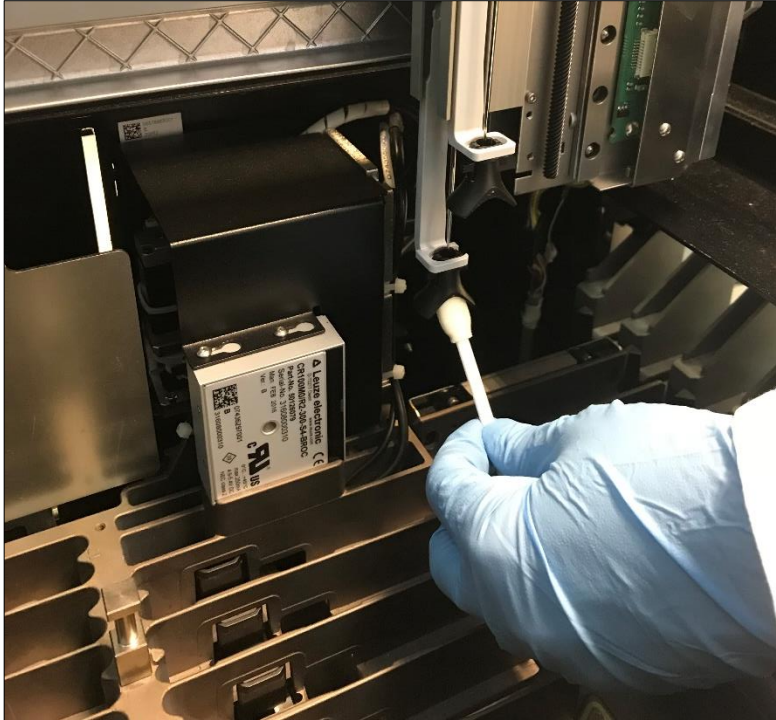

Now the maintenance tasks are completed, to finish the maintenance action, click on *Next*.

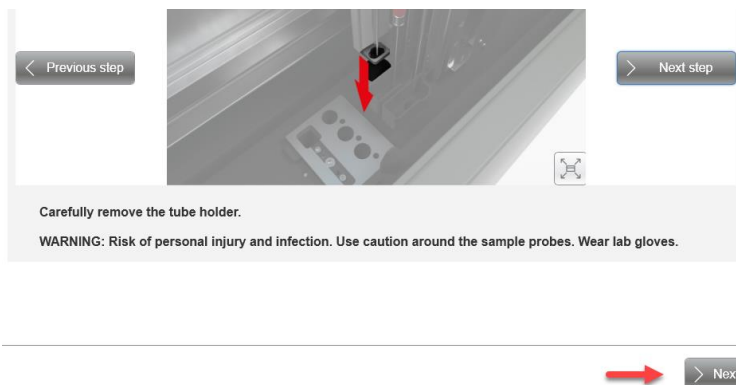

Follow the remaining instructions of the maintenance action.