

## Information Sheet

### VITROS® Chemistry Products AMYL Slides


#### VITROS® Chemistry Products Performance Verifier I and II Estimated SD

#### In Date Coatings of VITROS AMYL Slides:

Product Code (Unique Identifier)	Affected Coatings	GEN	Expiry Date (Range)
<b>1202670</b> (10758750008513) <b>8112724</b> (10758750004393)	0356	2	6/1/2020 – 9/1/2020
	0357	3	7/1/2020 – 2/1/2021
	0358	5	11/1/2020 – 7/1/2021
	0359	6	2/1/2021 – 10/1/2021
	0363	8	8/1/2021 – 10/1/2021
	0349	83	11/1/2020
	0362	84	5/1/2021
	0365	85	10/1/2021

#### Identifying Coating Numbers

Use the example below to determine the GEN and Coating ID:

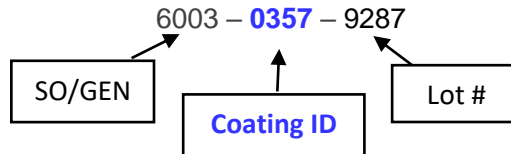

#### VITROS Performance Verifier I and II Estimated SD for VITROS AMYL Slides:

Level	GEN Range	Current Control Assay Sheet SD	Revised Control Assay Sheet SD	Units
PV I	1 - 49	4.2	6.8	U/L
	80 - 89	4.3	7.1	
PV II	1 - 49	10.5	17.1	
	80 - 89	10.8	17.6	